

Certificate of Analysis (Ver.2.0)

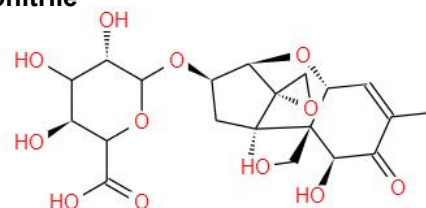
Deoxynivalenol-3-Glucuronide in Acetonitrile

1. General information

This document is designed and the certified value(s) and uncertainty(ies) are determined in accordance with ISO Guide 31^[1].

2. Description of the Reference Material (RM)

Product name:	Deoxynivalenol-3-Glucuronide in Acetonitrile
Product number:	STD#3107
CAS number:	1000000-13-4
Formula:	C ₂₁ H ₂₈ O ₁₂
Formula weight:	472.44
Lot#:	2B00D10
Result concentration:	25.06±0.46µg/mL
Starting material:	Deoxynivalenol-3-Glucuronide,lot#D21211P,Pribolab Pte. Ltd
Matrix:	Acetonitrile, LiChrosolv®, Merck
Amount:	1.2mL
Production date:	10,Mar,2022
Expiry date:	09,Sep,2023
Name of the supplier:	Pribolab Pte. Ltd



2.1 Intended use of the RM

- for laboratory use only
- calibration of analytical instruments

2.2 Instruction for the correct use of the RM

The compound should be stored at 2-8°C in a dark place. Before usage of the RM, the compound should be allowed to warm to temperature(20±3°C).The recommended minimum sub - sample amount for all kinds of application is 100 µL. The expiry date of this RM is based on the current knowledge and holds only for proper storage conditions in the originally closed flasks/ packages.

2.3 Hazardous situation

The normal laboratory safety precautions should be observed when working with this RM. Further details for the handling of this RM are available as safety data sheet.

Hazardous Ingredients	Concentration in%	Pictograms	Signal word	Hazard statement(s)
Acetonitrile	>99.9		Danger	H225,H302,H312,H319,H332

3. Certified values and their uncertainties

Deoxynivalenol-3-Glucuronide in Acetonitrile		
Compound	Mass concentration ^a	
Deoxynivalenol-3-Glucuronide	Certified value ^b	Uncertainty ^c
	25.06µg/mL	±0.46µg/mL
a Mass concentration based on weighed amount ,purity and dilution steps b Values are based on preparation data and confirmed experimentally by HPLC-DAD c Expanded uncertainty U(k=2) of the value u _c according to GUM ^[2]		

3.1 Calculation of uncertainty

After the concentration of the gravimetric prepared solution was confirmed by HPLC-DAD,the uncertainty of the calibrant was calculated on the basis of preparation^[3].

Uncertainty components	Description	Standard uncertainty(u)	
Purity(P)of solid Deoxynivalenol-3-Glucuronide	P=98.4±1.6%	u(P)=0.9%	a
Weighing procedure Weighted sample: m _{ws} =2.547mg	U _(m) =0.0000008g+1.30*10 ⁻⁵ *m _{Toxin} u _(m) =U _(m) /2	u(m)=0.0004mg	b
Dilution procedure Volumetric flask :V _f =100mL	Calibration:100mL ± 0.10mL Repeatability : 0.03mL Volume expansion solvent	u(cal)=0.04mL u(rep)=0.03mL u(Vol.exp.1)=0.24mL u(v)=0.25mL	c d e f
a Maximum tolerance of purity was divided by $\sqrt{3}$ b Calculation of this u-value is based upon the uncertainty formula for the weighed amount as given in the calibration report from annual balance calibration c A triangular distribution(division by $\sqrt{6}$)was chosen for the calculation of u(cal) d Based on a series of ten fill and weigh experiments on a typical 100mL flask; the value was used directly as a standard deviation e Based on the density of 0.7857 g/cm ³ at temperature T=20°C and a maximum temperature variation of ±3°C,of volume expansion, relative volume expansion coefficient of Acetonitrile is 1370*10 ⁻⁶ /°C[7],volume expansion term(rectangular distribution)was divided by $\sqrt{3}$ f The three contributions are combined to give the u(v)= $\sqrt{u(cal)^2 + u(rep)^2 + u(Vol.exp)^2}$			

Calculation of the combined uncertainty u_c and the expanded standard uncertainty U

$$C_{Toxin} = \frac{10 \times m_{ws} \times P}{V_f} = \frac{10 \times 2.547 \times 98.4}{100} = 25.06 \text{ mg / L}$$

$$\frac{u_c(C_{Toxin})}{C_{Toxin}} = \sqrt{\left[\frac{u(P)}{P}\right]^2 + \left[\frac{u(m)}{m_{ws}}\right]^2 + \left[\frac{u(V)}{V_f}\right]^2} = \sqrt{\left[\frac{0.9}{98.4}\right]^2 + \left[\frac{0.0004}{2.547}\right]^2 + \left[\frac{0.25}{100}\right]^2} = 0.009$$

$$u_c(C_{Toxin}) = C_{Toxin} \times 0.009 = 25.06 \times 0.009 = 0.23 \text{ mg / L}$$

calculation of expanded standard uncertainty U using a coverage factor k=2

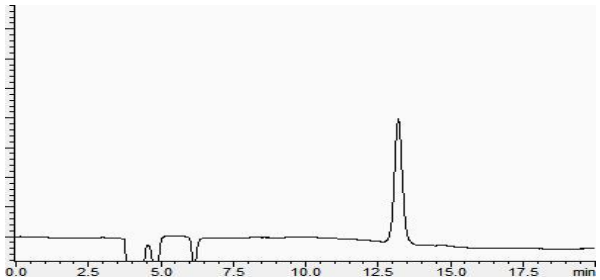
$$U(C_{Toxin}) = u_c(C_{Toxin}) \times 2 = 0.23 \times 2 = 0.46 \text{ µg / mL}$$

4. Discussion of traceability

This calibrant is certified on the basis of gravimetric preparation^[4]. Thus the certified value (mass concentration of Deoxynivalenol-3-Glucuronide) is based on the weighed amount of the starting material and is therefore traceable to the stated purity of the solid raw material. High purity material represents a practical realization of concentration units, through conversion of mass to molar quality.

5. Confirmation of certified value by HPLC-DAD

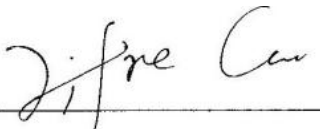
The certified concentration of Deoxynivalenol-3-Glucuronide of the gravimetric prepared solution was confirmed by HPLC-DAD against an independently prepared reference batch of Deoxynivalenol-3-Glucuronide.

column	T-3,150×4.6mm, 3.5µm	 <p>Figure 1: HPLC-DAD chromatogram of Deoxynivalenol-3-Glucuronide</p> <table border="1"> <thead> <tr> <th></th> <th>time [min]</th> <th>concentration ^a</th> </tr> </thead> <tbody> <tr> <td>Deoxynivalenol-3-Glucuronide</td> <td>13.200</td> <td>25.05 µg/mL</td> </tr> </tbody> </table> <p>^a Mean of 6 replicate measurements against reference batch, confidence interval with P=95%</p>		time [min]	concentration ^a	Deoxynivalenol-3-Glucuronide	13.200	25.05 µg/mL
	time [min]		concentration ^a					
Deoxynivalenol-3-Glucuronide	13.200		25.05 µg/mL					
injection Volume	20 µL							
solvent	acetonitrile / 0.1% formic acid in water = 10/90							
oven	30 °C							
flow rate	0.5 mL/min							
DAD settings	220 nm							
Sample dilution	acetonitrile/water = 10/90							

6. Further information

The purchaser must determine the suitability of this product for its particular use. Pribolab makes no warranty of any kind, express or implied, other than its products meet all quality control standards set by Pribolab. We do not guarantee that the product can be used for a special application.

Inspected by


 Quality System Specialist

References:

- [1] ISO Guide 31, 1-7, (2000), "Reference Materials - Contents of Certificates and Labels"
- [2] International Organization for Standardization (ISO), (2008), "Guide to the Expression of Uncertainty in Measurements ", (GUM 1995 with minor corrections) 1st Ed. Geneva, Switzerland
- [3] R.D. Josephs, R. Krska, S. MacDonald, P. Wilson, H. Pettersson, J. AOAC Int. 86, 50-60. (2003), "Preparation of a Calibrant as Certified Reference Material for Determination of the Fusarium Mycotoxin, Zearalenone"
- [4] E.W. Flick, (1998), "Industrial Solvents Handbook ",5rd Ed., Noyes Data Corp. Westwood NJ