# Certificate of Analysis(Ver.1.0)

U-[<sup>13</sup>C<sub>20</sub>]-T-2 triol in Acetonitrile

## 1. General information

This document is designed and the certified value(s) and uncertainty(ies) are determined in accordance with ISO Guide 31[1].

# 2. Description of the Reference Material (RM)

Name:	U-[ <sup>13</sup> C <sub>20</sub> ]-T-2 triol in Acetonitrile
Catalog number:	STD#3112U
CAS number:	34114-98-2(unlabeled)
Formula:	<sup>13</sup> C <sub>20</sub> H <sub>30</sub> O <sub>7</sub>
Formula weight:	402.45
Lot #:	1101A23
Starting material:	U-[ <sup>13</sup> C <sub>20</sub> ]-T-2 triol,lot#J19117P,Pribolab Pte.Ltd.
Solvent:	Acetonitrile,LiChrosolv <sup>®</sup> ,Merck
Amount:	1.2mL
Production date:	25/10/2020
Expiry date:	24/10/2021
Name of the supplier:	Pribolab Pte.Ltd.

# 2.1 Intended use of the RM

- for laboratory use only

- internal standard[2]

#### 2.2 Instruction for the correct use of the RM

The compound should be stored at 2-8°C in a dark place. Before usage of the RM , the compound should be allowed to warm to temperature( $20\pm3$ °C). The recommended minimum sub-sample amount for all kinds of application is 100 µL . The expiry date of this RM is based on the current knowledge and holds only for proper storage conditions in the originally closed flasks/packages.

# 2.3 Hazardous situation

The normal laboratory safety precautions should be observed when working with this RM.Further details for the handing of this RM are available as safety data sheet.

Hazardous Ingredients Acetonitrile Concentration in% >99.9



Signal word Danger Hazard statement(s) H225,H302,H312,H319,H332

Pribolab Pte.Ltd.Singapore. 10 Biopolis Road #15-06 Immunos 138670 Singapore www.pribolab.com

#### 3. Certified values and their uncertainties

U-[ <sup>13</sup> C <sub>20</sub> ]-T-2 triol in	acetonitrile	
Compound	Mass conce	entration <sup>a</sup>
U-[ <sup>13</sup> C <sub>20</sub> ]-T-2 triol ,98.6 atom% <sup>13</sup> C	Certified value <sup>b</sup>	Uncertainty <sup>c</sup>
0-[**C <sub>20</sub> ]-1-2 (101,98.0 a(011% **C	10.03µg/mL	±0.14µg/mL
a Values are based on preparation data and confi	rmed experimentally by H	IPLC-DAD
b Mass concentration based on weighed amount,	purity and dilution step	
c Expanded uncertainty $U(k=2)$ of the value $u_c$ acc	ording to GUM[3]	

# 3.1 Calculation of uncertainty

The uncertainty of the calibrant solution was calculated on the basis of preparation [4].

Uncertainty components	Description	Standard uncertainty(U)	
Purity(P)of solid	P=98.8±1.2%	u(P)=0.7%	а
U-[ <sup>13</sup> C <sub>20</sub> ]-T-2 triol ,98.6 atom% <sup>13</sup> C			ľ
Weighing procedure	U <sub>(m)</sub> =0.0000008g+1.30*10 <sup>-5</sup> *m <sub>Toxin</sub>	u <sub>(m)</sub> =0.0004mg	b
Weighted sample: m <sub>ws</sub> =1.015mg	u <sub>(m)</sub> =U <sub>(m)</sub> /2	u(m)=0.0004mg	
	calibration: 100±0.1mL	u(cal)=0.04mL	С
Dilution procedure	repeatability: 0.04mL	u(rep)=0.04mL	d
Volumetric flask:V <sub>f</sub> =100mL	volume expansion solvent	u(Vol.exp.)=0.24mL	e
		u(v)=0.2mL	f

a Maximum tolerance of purity (rectangular distribution) was divided by  $\sqrt{3}$ 

b Calculation of this u-value is based upon the uncertainty formula for the weighed amount as given in the calibration report from annual balance calibration

c A triangular distribution ( division by  $\sqrt{6}$  ) was chosen for the calculation of u(cal)

d Based on a series of ten fill and weigh experiments on a typical 100mL flask; the value was used directly as a standard deviation

e Based on the density of 0.7857 g/cm<sup>3</sup> at temperature T=20°C and a maximum temperature variation of  $\pm$ 3°C, of volume expansion, relative volume expansion coefficient of acetonitrile is 1370\*10<sup>-6</sup>/°C[5], volume expansion term(rectangular distribution) was divided by  $\sqrt{3}$ 

f The three contributions are combined to give the  $u(V) = \sqrt{u(cal)^2 + u(rep)^2 + u(Vol.exp)^2}$ 

Calculation of the combined uncertainty u c and the expanded standard uncertainty U

$$C_{\text{Toxin}} = \frac{10 \times \underline{m}_{\text{ws}} \times P}{V_{f}} = \frac{10 \times 1.015 \times 98.8}{100} = 10.03 \text{ mg/L}$$
$$\frac{u_{c}(c_{\text{Toxin}})}{c_{\text{Toxin}}} = \sqrt{\left[\frac{u(P)}{P}\right]^{2} + \left[\frac{u(m)}{m_{\text{WS}}}\right]^{2} + \left[\frac{u(V)}{V_{f}}\right]^{2}} = \sqrt{\left[\frac{0.7}{98.8}\right]^{2} + \left[\frac{0.0004}{1.015}\right]^{2} + \left[\frac{0.2}{100}\right]^{2}} = 0.007$$

$$\mathbf{u}_{\rm C} (\mathbf{C}_{Toxin}) = \mathbf{C}_{Toxin} \times 0.007 = 10.03 \times 0.007 = 0.07 \, mg \, / L$$

calculation of expanded standard uncertainty U using a coverage factor k=2

$$U(\mathbf{C}_{Toxin}) = \boldsymbol{\mu}_{c}(\boldsymbol{C}_{Toxin}) \times 2 = 0.07 \times 2 = 0.14 \, \mu g \, / \, mh$$

Pribolab Pte.Ltd.Singapore. 10 Biopolis Road #15-06 Immunos 138670 Singapore www.pribolab.com

Isotope pat	tern <sup>a</sup>
Compound	Isotopic distribution
U-[ <sup>13</sup> C <sub>20</sub> ]-T-2 triol	77.8%
U-[ <sup>13</sup> C <sub>19</sub> ]-T-2 triol	15.4%
U-[ <sup>13</sup> C <sub>18</sub> ]-T-2 triol	6.8%
Calculated isotopic enrichmer	nt levelª:98.6atom% <sup>13</sup> c
<sup>a</sup> Approximation based on LC-MS/MS data	


## 4.Isotopic enrichment and isotope pattern

#### 5. Discussion of traceability

This calibrant is certified on the basis of gravimetric preparation [5]. Thus the certified value(mass concentration of U-[ $^{13}C_{20}$ ]-T-2 triol ,98.6atom%  $^{13}C$  is based on the weighed amount of the starting material and is therefore traceable to the stated purity of the solid raw material. High purity material represents a practical realization of concentration units, through conversion of mass to molar quantity.

# 6. Confirmation of certified value by HPLC-DAD

The certified concentration of U-  $[^{13}C_{20}]$ -T-2 triol,98.6 atom%  $^{13}C$  of the gravimetric prepared solution was confirmed by HPLC-DAD against an independently prepared reference batch of unlabeled T-2 triol.

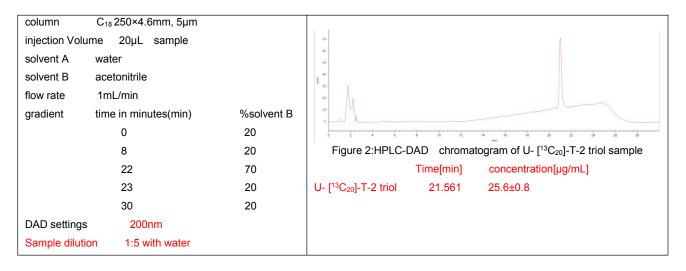

Pribolab Pte.Ltd.Singapore. 10 Biopolis Road #15-06 Immunos 138670 Singapore www.pribolab.com

#### 7. Further information

The purchaser must determine the suitability of this product for its particular use. Pribolab makes no warranty of any kind, express or implied, other than its products meet all quality control standards set by Pribolab. We do not guarantee that the product can be used for a special application.

Inspected by

re

Quality System Specialist

#### **References:**

- [1]ISO Guide 31:2015 1-18, "Reference materials contents of certificates, labels and accompanying documentation"
- [2]G. Häubl, F. Berthiller, R. Krska, R. Schuhmacher, "Suitability of a fully <sup>13</sup>C isotope labelled internal standard for the determination of the mycotoxin deoxynivalenol by LC-MS/MS without clean-up", Anal. Bioanal. Chem. 384 (3), (2006), 692-696
- [3] International Organization for Standardization (ISO), (2008), "Guide to the expression of uncertainty in measurement", (GUM 1995 with minor corrections) 1<sup>st</sup> Ed. Geneva, Switzerland
- [4] R.D. Josephs, R. Krska, S. MacDonald, P. Wilson, H. Pettersson, J. AOAC Int. 86, 50-60, (2003),
  "Preparation of a Calibrant as Certified Reference Material for Determination of the Fusarium Mycotoxin Zearalenone"
- [5] E.W. Flick, (1998), "Industrial Solvents Handbook", 5th Ed., Noyes Data Corp. Westwood NJ